

# BioPlus 1270/1400



Biostorage you can depend on



BioPlus 1400.  
Drawers and glass doors are optional extras.


The **BioPlus 1270/1400** range is designed for the storage of the most delicate biomaterials, in situations where even tiny fluctuations in conditions inside the storage cabinet can have a serious effect on the contents.

BioPlus cabinets also enable you to reduce internal relative humidity, thus cutting down on any likelihood of undesirable contaminants getting near delicate biomaterials.

In terms of performance, these units are designed to provide the very best results even under exceptional conditions. The BioPlus utilises environmentally responsible HFC-free foam propellants and is available with HFC-free refrigerants.



The BioPlus 1270/1400 is available as a refrigerator or freezer in either white or stainless steel finish. All BioPlus models come with self-closing door. Glass door is available for refrigerator.

 II 3G Ex nA nC nL IIB T3 Gc

Biostorage  
you can  
depend on

Version 06 Valid from 01.01.23



# BioPlus 1270/1400



Biostorage you can depend on

## Top panel with BioLine control unit

The MPC unit is specially developed for users of BioLine storage cabinets. Among its many features are acoustic temperature and door alarms, and a voltage-free contact (NO/NC) for remote alarms. The display can be easily locked.




## Access port

24.5 mm port for ease of access, e.g. sensors for external temperature surveillance.



## E-sensor

Extra sensor for providing a temperature reference within the storage compartment.

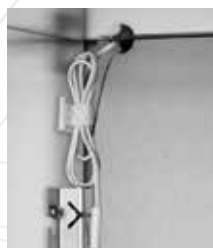

## Shelves and drawers

A wide range of versatile interior fittings and layouts is available. The operating specifications and the interior layout can both be customised.




## Stainless steel interior

The internal sides and bottom as well as shelf and drawer brackets and support rails are made of stainless steel.

## ATEX

All BioPlus models comply with ATEX Cat.3 Zone 2 requirements, both internal and externally, allowing for placement in locations where explosive gas atmospheres may be present.

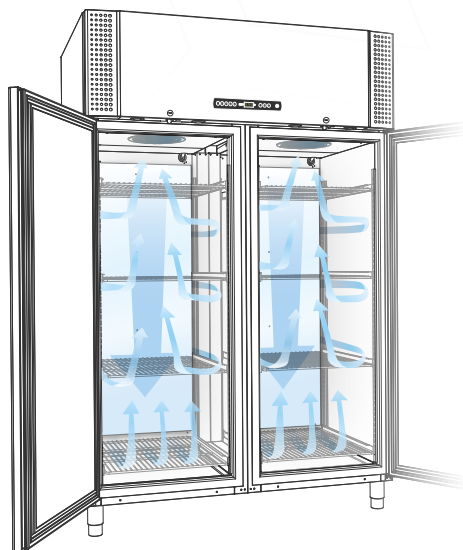

## AIR DISTRIBUTION

The unique Gram BioLine air distribution system makes sure the temperature inside the cabinet remains stable at all times, while using the minimum amount of energy.

Technical specifications	BioPlus ER1270/1400	BioPlus RF1270/1400
Temperature ranges	+20/-2 °C	-5/-25 °C
Ambient temperature range	Solid door +10/+43 °C; glass door (ER) +10/+38 °C	
Control Unit	Gram Control Unit with voltage-free contact, E-sensor, acoustic and visual door and temperature alarms that can be programmed individually, alarm recording and offset function	
Material interior	Stainless steel	
Material exterior	White lacquered steel or stainless steel	
Dimensions mm (W x D x H)	1270: 1390 x 876 x 1875/2125 1400: 1390 x 876 x 2025/2275	
Gross volume	1270: 1270 litres/44.8 cubic feet / 1400: 1400 litres/49.4 cubic feet	
Net volume	1270: 864 litres/30.6 cubic feet / 1400: 968 litres/34 cubic feet	
Net weight	Unpacked net weight without optional fittings: 1270: 271 kg / 1400: 283 kg	
Modules for shelves and drawers	1270: Shelves 48 - Drawers 24 1400: Shelves 54 - Drawers 26	
Insulation	60 mm polyurethane with HFC-free cyclopentane propellant	
K-Value	Glass door 0.50/0.50 W/m²K	Solid door 0.31/0.31 W/m²K
Refrigerant	R290	
CO <sub>2</sub> e	-	
Refrigeration capacity at -10 °C	R290 946 watt	
Refrigeration capacity at -25 °C	R290 950 watt	
Energy consumption	2.50 kWh/24h/ 2.56 kWh/24h	7.97kWh/24h/ 8.39 kWh/24h
Base	Legs**	
Sound level	48.3 dB(A) / 47.5 dB(A)	51 dB(A) / 46.6 dB(A)
Connection	230 V, 50 Hz	
Air system	Gram BioLine ventilated air distribution system	
Defrost system	Automatic smart defrost with re-evaporation of defrost water	

\*\* The cabinets can also be delivered with castors or for mounting on a plinth at no extra charge